

RNAscope[®] Chromogenic Assay for PBMC and Non-Adherent Cells (Part 1)

Introduction

This Technical Note provides guidelines for the preparation of Peripheral Blood Mononuclear Cells (PBMC) and non-adherent cells that can be assayed using an RNAscope[®] Chromogenic Detection Kit. The required RNAscope[®] Pretreat Reagents are Pretreat 1 (Cat. No. 310020) and Pretreat 3 (available in Cat. No. 320842). In addition, the following reagents are required RNAscope[®] PBMC Preparation Reagents (Cat. No. 320970, includes Cell Prep and PBMC Wash

reagents). In addition material required but not provided by ACD–100% EtOH, Histopaque 1077 (Sigma-Aldrich), SuperFrost plus microscope slide (Fisher) and PBMC Fix (10%NBF). For every chemical, read the Material Safety Data Sheet (MSDS) and follow handling instructions. For every chemical, read the Material Safety Data Sheet (MSDS) and follow handling instructions.

For the latest service and support information, go to: www.acdbio.com/support

Part 1: Prepare Samples

PBMC Collection

Reagent Preparation


1. Add 3 mL of Histopaque[®] 1077 (or other Ficoll solution) to a 15 mL conical centrifuge tube and bring the solution to **ROOM TEMPERATURE (RT)**.
2. Prepare CellPrep/PBS solution by dissolving 50X CellPrep (stock should be frozen at **-20°C**) in 1X PBS. Prepare 50 mL for each 5 mL blood sample.

PBMC Purification

1. Transfer 5 mL of blood to a 15 mL conical tube.
2. Carefully layer the 5 mL of blood sample onto the Histopaque[®] 1077 solution.
3. Centrifuge at **RT** in a horizontal rotor (swing-out head) for **20 MIN** at 800 RCF (with minimal acceleration/break).
4. Carefully remove the upper phase (plasma phase) with a pipette or aspiration device, leaving ~0.5 cm above the PBMC layer.
5. Use a 1 mL pipette to transfer the PBMC layer to the 50 mL polypropylene tube containing 40 mL

CellPrep/PBS. Pipette up and down several times to minimize cell loss in the pipette tip.

6. Centrifuge at **RT** for **10 MIN** at 250 RCF (with maximum acceleration/break).
7. Aspirate supernatant without disturbing the cell pellet, leaving ~ 5 mL liquid.
8. Resuspend cell pellet with remaining liquid by pipetting up and down 10 times then transfer to a

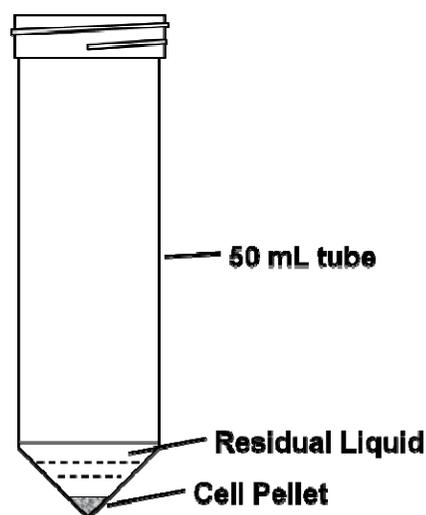



new 15 mL tube.

9. Wash the 50 mL tube with 5 mL CellPrep/PBS solution and transfer the solution to the 15 mL tube containing resuspended cells to minimize cell loss.
10. Centrifuge at **RT** at 250 RCF for **10 MIN** (with maximum acceleration/break).
11. Aspirate supernatant leaving as little liquid as possible without touching the cell pellet.

Non-adherent Cell Collection

1. Harvest cells by centrifuging at **RT** at 250 RCF for **10 MIN** in a 50 mL polypropylene tube.
2. Aspirate supernatant without disturbing the cell pellet.
3. Wash with 40 mL 1X PBS by resuspending cells and centrifuging at RT at 250 RCF.
4. Aspirate supernatant, leaving as little liquid as possible without touching the cell pellet

Cell Fixation

1. Resuspend cells in 5 mL of PBMC-Fix. Gently pipette up and down 10 times to completely break apart the cell pellet.
2. Incubate the tube in a **37°C** water bath for **1 HR**.

Post-Fixation Wash and Storage

1. Centrifuge at 250 RCF for **10 MIN** to pellet the cells.
2. Remove supernatant without disturbing the pellet.
3. Resuspend cells in 10 mL PBMC-Wash, and centrifuge at 250 RCF for **10 MIN**.
4. Resuspend cells in 10 mL PBMC-Preserve. Pipette up and down 10 times to completely break apart the cell pellet.
5. Incubate at **RT** for **10 MIN** and transfer to **4°C**.

NOTE: The cells can be stored in PBMC-Preserve at **4°C** for up to **7 days**.

Slide Preparation

1. Adjust the cell density with PBMC-Preserve (70% EtOH) to 1×10^6 cells/mL.

2. Mix well by pipetting. Transfer 1 mL cell suspension to each pre-assembled cyto-centrifuge cartridge.

NOTE: Cell density and volume described here is based on the Hettich cyto-centrifuge with an 8 mL funnel chamber. If other cyto-centrifuge systems are used, adjust the cell density and volume to achieve a single cell layer after cyto-centrifuge.

3. Cyto-centrifuge at 800 RCF for **20 MIN**.
4. Carefully remove supernatant completely with pipette, disassemble cyto-centrifuge cartridge.
5. Air dry slides for **20 MIN** at **RT**.
6. Immerse slides in 50% EtOH. Incubate at **RT** for **5 MIN**.
7. Remove 50% EtOH and replace with 70% EtOH. Incubate at **RT** for **5 MIN**.
8. Remove 70% EtOH and replace with 100% EtOH. Incubate at **RT** for **5 MIN**.
9. Remove 100% EtOH and replace with fresh 100% EtOH. Incubate at **RT** for **5 MIN**.

NOTE: The slides can be stored in 100% EtOH at **-20°C** for up to **1 MONTH**.

Part 2: RNAscope® Pretreatment

Prepare Materials

1. Bring HybEZ™ Oven to **40°C**.
2. Place a wet humidifying paper in the Humidity Control Tray, leaving the HybEZ™ Slide Rack on bench. Re-insert the covered tray into the oven and close the oven door. The tray should be pre-warmed for at least **20 MIN** before use.



Create a Hydrophobic Barrier

1. Remove slides from 100% EtOH and dry at **37°C** for **30 MIN** on a slide warmer.
2. Draw 2–4 times around the cell spot using the Immedge™ hydrophobic barrier pen. Let the barrier dry completely **~1 MIN**.

Add Pretreat 1

1. Add 2–4 drops Pretreat 1 for **10 MIN** at **RT**. Use enough solution to completely cover the cell spot.

NOTE: Avoid dropping solution directly onto the cell spot to minimize cell loss. Instead, add reagent around the edge of the cell spot.

2. One at a time, tap/flick each slide to remove excess liquid and submerge in 1X PBS.
3. Remove 1X PBS, replace with fresh 1X PBS, and incubate at **RT** for **1 MIN**.

Add Pretreat 3

1. One at a time, remove each slide from the 1X PBS and tap/flick to remove excess liquid. Place the slides on the HybEZ™ Slide Rack.
2. Add 2–4 drops Pretreat 3. Use enough solution to completely cover the cell spot.

For Molecular Biology Applications (MBA), not intended for diagnosis. Refer to appropriate regulations.

NOTICE TO PURCHASER: PLEASE REFER TO THE RNASCOPE® 2.0 ASSAY USER MANUAL FOR LIMITED USE LABEL LICENSE OR DISCLAIMER INFORMATION. Advanced Cell Diagnostics, Inc. reserves the right to change its products and services at any time to incorporate technological developments. This manual is subject to change without notice. Although this manual has been prepared with every precaution to ensure accuracy, Advanced Cell Diagnostics, Inc. assumes no liability for any errors, omissions, or for any damages resulting from the use of this information.

© 2014 Advanced Cell Diagnostics. All rights reserved. RNAscope® and HybEZ™ are trademarks of Advanced Cell Diagnostics, Inc. All other trademarks belong to their respective owners.

3. Place the slide rack in the pre-warmed Humidity Control Tray, close lid, and incubate the tray in the HybEZ™ Oven for **30 MIN** at **40°C**.
4. Take slides out of the oven and one at a time tap/flick to remove excess Pretreat 3. Do not let sample dry out.
5. Submerge the slides in a Coplin jar containing 1X PBS.
6. Remove 1X PBS, replace with fresh 1X PBS, and incubate at **RT** for **1 MIN**.

IMPORTANT! Proceed to the RNAscope® protocol using the appropriate Part 2 Detection User Manual* available at www.acdbio.com/support/technical-doc.

* RNAscope® 2.0 HD Detection Kit-Brown User Manual, Part2 (Cat. No.320497); RNAscope® 2.0 HD Detection Kit-Red User Manual, Part 2 (Cat. No. 320487); RNAscope 2-Plex Detection Kit-Chromogenic User Manual (Cat. No.320494)

Obtaining Support

For the latest services and support information, go to: www.acdbio.com/support.

At the website, you can:

- Access telephone and fax numbers to contact Technical Support and Sales.
- Search through FAQs.
- Submit a question directly to Technical Support.

Unlock your RNA with RNAscope® Technology

RNAscope® allows single-cell gene expression profiling *in situ*—unlocking the full potential of RNA biomarkers. The targeted molecular signature of every cell in a sample is revealed and measured precisely, all within the intricate cellular and tissue architecture of clinical specimens.